



Nucleic Acid Rapid Test Vehicle

Efficient and Safe Test, Ready at Any Time



- ✓Antigen Sampling + Test
- ✓Nucleic Acid Sampling + Test

- **Special Designed:** The whole vehicle is special designed and modified, with reasonable layout, and conforms with national vehicle standards.
- **Professional and Ensured:** Biological safety cabinet ensures independent space for samples and medical personnel to avoid cross-infection
- **Multi-function:** It supports antigen test, nucleic acid sampling and rapid test, flexible and changeable.
- **Safety:** Positive and negative pressure control system, H14 level filtration and exhaust system, ultraviolet disinfection system
- **High Test Efficiency:** High Test efficiency, the results are available within a short period of time
- **Intelligent Control:** Automatic voltage stabilization, abnormal condition alarm, automatic disinfecting



Qingdao Haier Biomedical Co.,Ltd.

No.280 Feng Yuan Road, High-tech Zone,
Qingdao, 266109, P.R. China
Tel: +86-0532-88935593
E-mail: inquiry@haierbiomedical.com
Website: www.haiermedical.com



Haier Biomedical
International



Haier Biomedical
International



@haiermedicalint



Haier Biomedical
International



Haier Biomedical
International

Note: If a slight difference occurs between pictures and actual products, please refer to actual products. Our company reserves the right of final interpretation of this brochure, please contact us for any further information if required.

*Haier Biomedical reserves the right to change products and specifications without prior notice.



Nucleic Acid Rapid Test Vehicle

Usage/Scope of Application

- Suitable for towns and rural areas where people are scattered; and communities, expressway intersections, airports, public places and assemblies where emergency sampling and rapid test for nucleic acid and antigen test are temporarily required.

Product Advantages

➤ High-efficiency Air Filtration System

- H14 filtration and exhaust system, the filtration efficiency of 0.3μm particles reaches 99.999%, ensuring clean exhaust air
- Intelligent ultraviolet disinfecting system, delayed opening and timed closing, safer disinfecting operation

➤ Efficient Disinfection System


- Equipped with UV lamp which has delayed opening and timed closing function
- Infrared induction for personnel protection

➤ Comfortable and Safe Environment

- IoT safety cabinet, negative pressure in the test area, ensuring the safety of materials, people and the environment
- Fully enclosed internal and external interlocking transfer window, built-in UV disinfection system, ensuring the safety of sample transfer

➤ Other Product Functions

- Energy-saving heating and cooling air conditioning, remote temperature adjustment, more comfortable working environment
- Vehicle-mounted monitoring station, high mobility, ready at any time

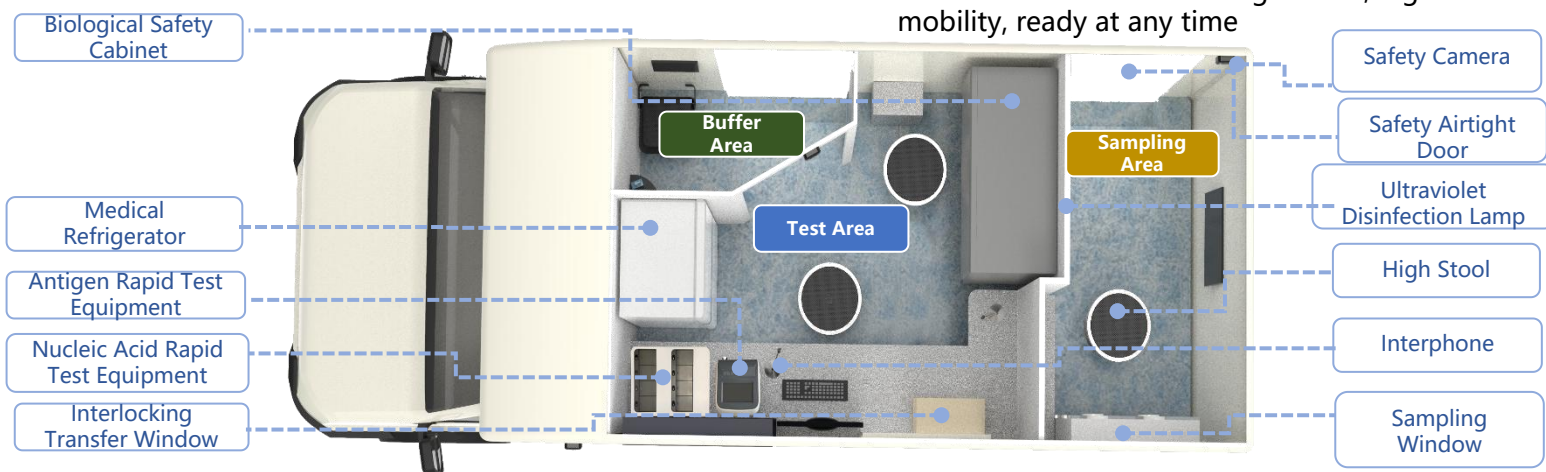

Product Specifications

L	W	H	Seat	Max. Speed	Fuel Consumption	Fuel Tank
5995mm	2200mm	2865mm	3 seats	145KM/H	10L/100 km	100L
Vehicle Weight	Drive Mode	Rated Power	Filtration Level	Chassis		Engine Displacement
4250KG	4x2 front rear drive	95kW	H14	IVECO Daily Series		2.287L

Qingdao Haier Biomedical Co.,Ltd.

No.280 Feng Yuan Road, High-tech Zone,
Qingdao, 266109, P.R. China
Tel: +86-0532-88935593
E-mail: inquiry@haierbiomedical.com
Website: www.haiermedical.com



Haier Biomedical
International



Haier Biomedical
International



@haiermedicalint



Haier Biomedical
International



Haier Biomedical
International

Note: If a slight difference occurs between pictures and actual products, please refer to actual products. Our company reserves the right of final interpretation of this brochure, please contact us for any further information if required.

*Haier Biomedical reserves the right to change products and specifications without prior notice.